



*Tablet and Capsule  
Counting Solutions*



## Technology You Can Rely On

The electronic tablet and capsule counters in the Romaco Bosspak RTC series impress with their unique optical matrix counter. Quad Count ensures extremely precise filling processes by scanning every single product from four sides. Broken tablets and other defective products are reliably identified. The technology is therefore particularly recommended for processing high quality and potent pharmaceutical solids. With a maximum counting speed of 50 bottles a minute, Bosspak counters are mainly designed for small to medium batch sizes.

## Quad Count Optical Matrix Counter

- Unique counting system scans products from four sides
- Larger count zone provides very high counting accuracy
- Fail-safe monitoring of the counting sensors
- Special nozzles prevent dust on the counting heads
- Broken and defective products are reliably detected
- Suitable for validating conventional counters

## Reliable Filling

- Products separated properly by means of rotating sorting discs
- Single-stream filling (products fall one by one into the bottle)
- Vibration-free filling
- Efficient dust aspiration before counting and filling
- Special coating prevents static on the products

## User Friendly Design

- Quick, tool-free format changes
- Easily accessible and simple to clean
- Low noise level
- Fully or semi-automatic operation


## Versatile Applications

- Handling of tablets, soft or hard gelatine capsules, dragées and caplets in a variety of shapes and sizes
- Even very small quantities can be counted (from one product per bottle)
- cGMP compliant design
- Suitable for small to medium batch sizes

*Exploiting Technologies  
Romaco Bosspak RTC Series*



*Product separation*

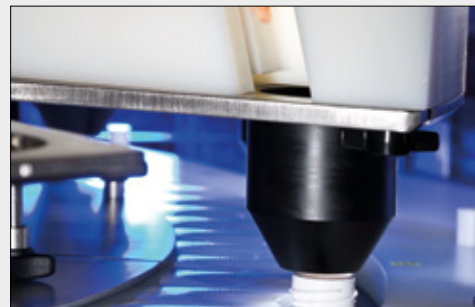

*Bottling*



*Primary packaging for pharmaceutical solids*



**Romaco Pharmatechnik GmbH**

Am Heegwald 11  
76227 Karlsruhe  
Germany

P +49 (0)721 4804 0  
F +49 (0)721 4804 211  
E [bosspak@romaco.com](mailto:bosspak@romaco.com)

[www.romaco.com](http://www.romaco.com)

 **INNOJET**  
*Herbert Hüttlin*

 **KILIAN**

 **NOACK**

 **SIEBLER**

 **PROMATIC**

 **MACOFAR**

 **BOSSPAK**

