


Bottles Depacking system MED 003 from Romaco Macofar

The Romaco Macofar MED 003 is a depacking machine for bottles packed in thermos-shrunk plastic. The machine automatically cuts and remove the thermos-shrunk plastic, freeing the bottles from their package. Thanks to the brushless motorization, all operations are smooth and maintenance needs are very low.

Highlights

- Compact footprint
- Brushless motorization granting repeatability of process
- Receipts recording
- Tool-less changeover
- Suitable for full automation of loading operations
- Excellent overview of, and accessibility to, the entire system for easy maintenance and setting according to Poka-Yoke philosophy
- Automatic ejection of plastic film into dedicated rejection bin
- Thermos-shrunk plastic material cut on 3 sides by means of thermal blade

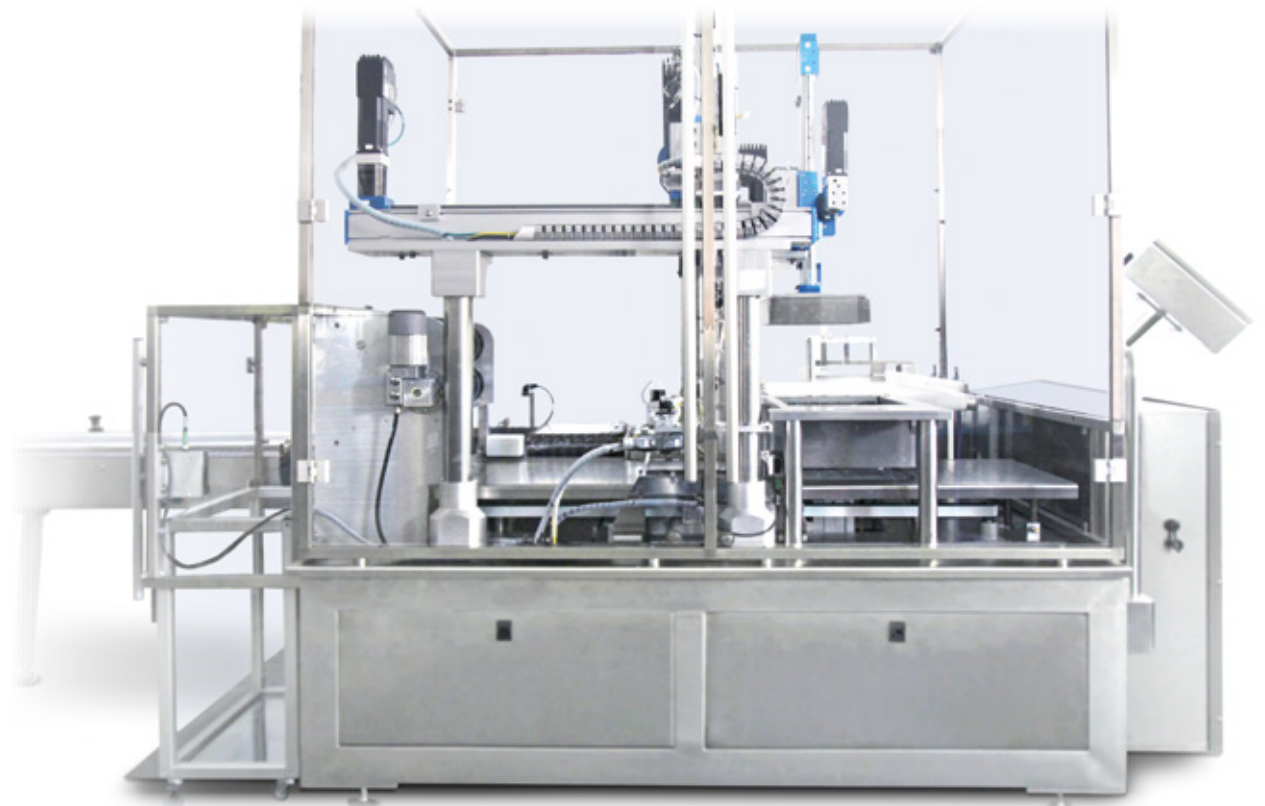

Technical Data	MED 003	
Mechanical speed (cpm) approx.	2.5	
Weight (kg) approx.	700	
Length (mm)	2,020 plus inlet belt (total 3,460)	
Depth (mm)	1,785 plus outlet belt (total 2,700)	
Height (mm)	1,950	
Vials size range (mm)	min.	max.
Diam.	20	100
H.	40	200
Bundle size range (mm)	min.	max.
Width	200	500
Length	300	600
Height	40	200

