


## HM 1E series from Romaco Siebler

The Siebler HM 1E series is specifically designed for heat sealing effervescent tablets, a highly sensitive product that requires the best possible protection against humidity. The sealing parameters comprising pressure, temperature and speed are pre-set and stored as format parameters for the relevant products. For contract manufacturing companies or customers with different pack sizes the new sealing system with pluggable tubular format tools is a perfect solution.

### Highlights

- High availability
- Quick production start
- Top sealing quality (repeatable/pre-set)
- Reproducible sealing parameters
- Perfect accessibility
- Ready for Remote Assist
- OPC UA interface for production data acquisition (PDA)

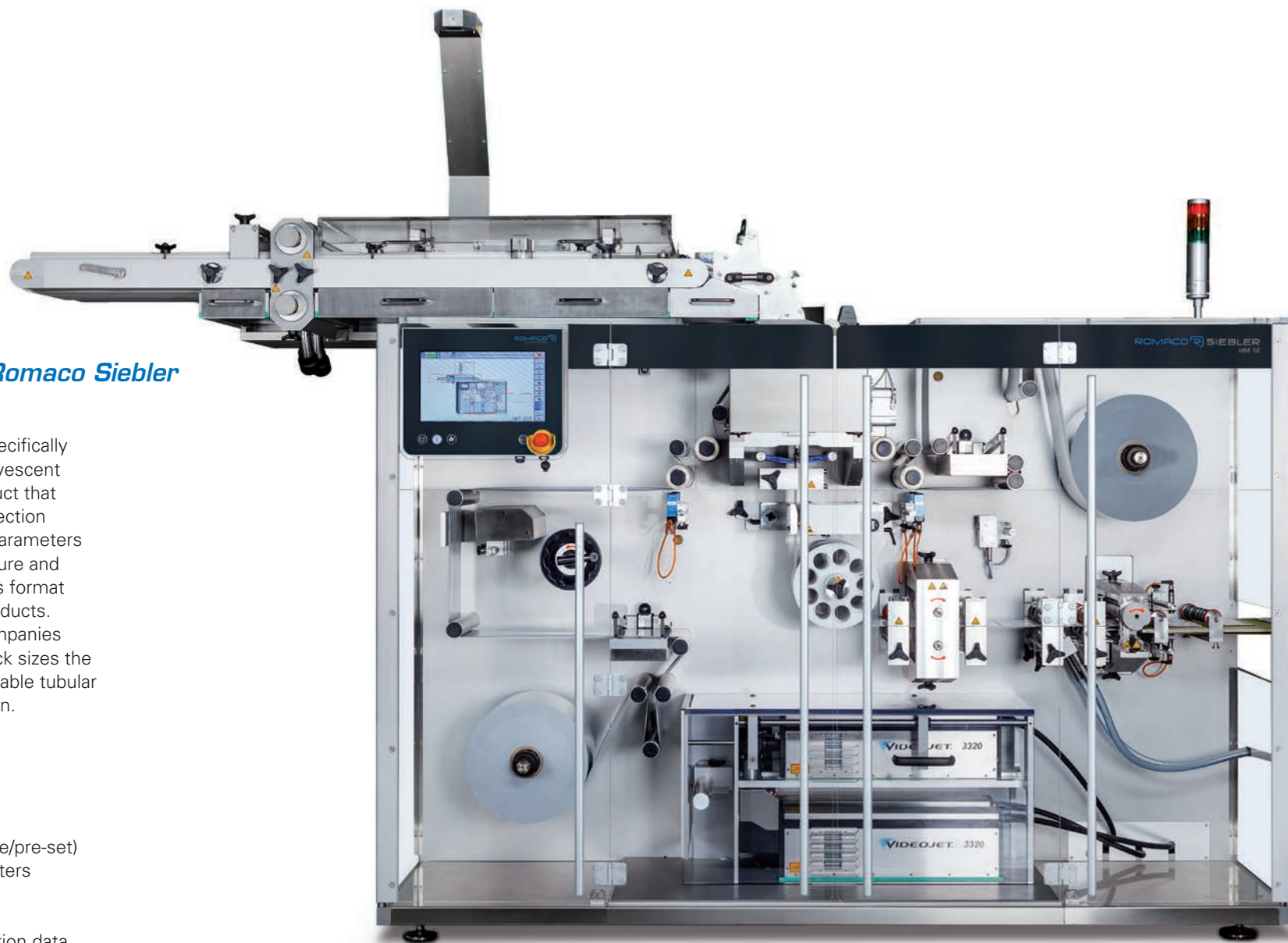

<b>Technical Data</b>	HM 1E-350	HM 1E-600
Foil width max., in mm	360	600
Product lanes max.	8	12
Output tablets/minute max.	2,000	3,000
Remote Assist connection	Ethernet	
PDA interface	OPC UA	

Technical changes reserved

