

Romaco Siebler HM 3


The Siebler HM 3 heat-sealing machines produce sachet packs for pharmaceutical granules. The air, light and moisture-tight, four-sided sealed foil packs are manufactured at very high speeds. Line configurations with Innojet processing machines, Siebler transfer solutions and Promatic secondary and tertiary packaging machines can be implemented at any time.

Highlights

- GMP compliant design with balcony architecture
- Very compact machine with a small footprint
- Feeding unit also optimally accessible
- Excellent OEE values
- Remote Assist with Call 4 Service, smart glasses and status app

Benefits

- Premium sealing quality thanks to QuickSeal
- Servo driven product feeding
- Highly accurate dosing
- Multiple dosing of two granules into one sachet
- Double packs with two different products
- Print mark registration for both foil rolls

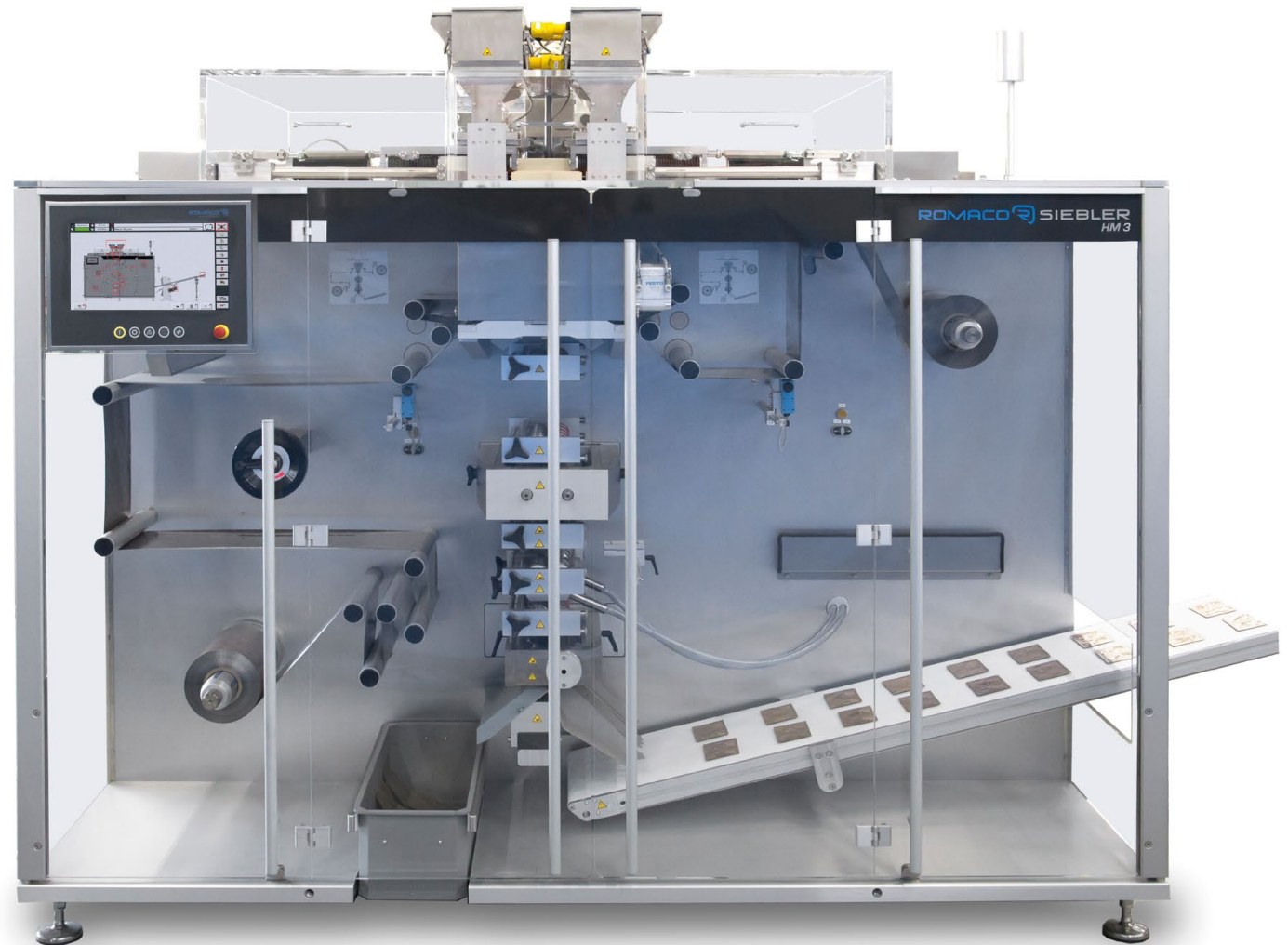

Technical Data	HM 3-350
Foil width max., in mm	360
Product lanes max.	8
Output sachets/minute max.	960
Remote Assist connection	Ethernet
PDA interface	OPC UA

Technical changes reserved
Other designs on request

